


Sample Information

Name : 二肽-6
Sequence : Pro-Hyp
Pump A : 0.1%Trifluoroacetic in 100% water
Pump B : 0.1%Trifluoroacetic in 100% acetonitrile
Total Flow : 1ml/min
Wavelength : 214nm
Analytical column type : Phenomenex Luna® 5 µm C18(2) (4.6*250mm*5 µm)
Dissolution method : 100%H₂O
Inj. Volume : 25ul

Time	Module	Action	Value
0.01	Pumps	B.Conc	0.1
20.00	Pumps	B.Conc	20
33.00	Pumps	B.Conc	100
38.00	Pumps	B.Conc	100
40.00	Pumps	B.Conc	0.1
50.00	Controller	Stop	

Chromatogram

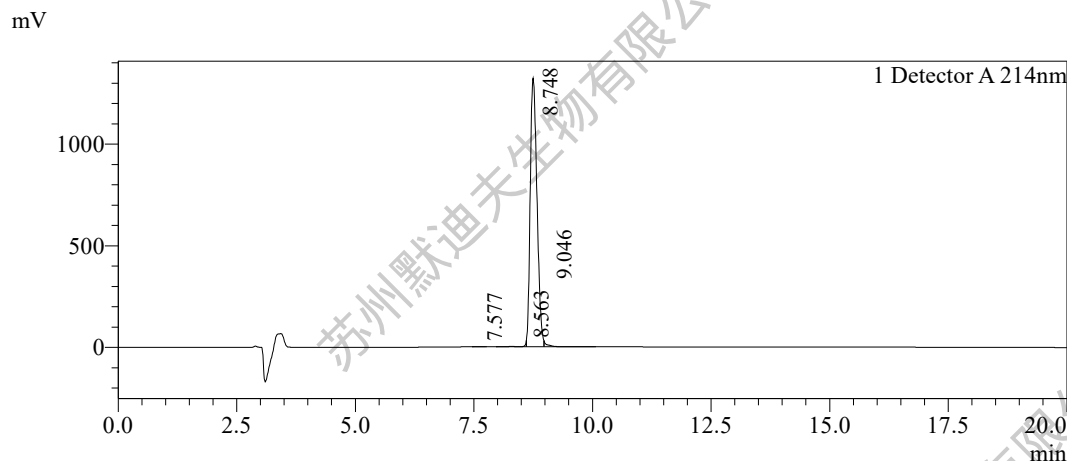

Peak Table

Detector A 214nm

Peak#	Ret. Time	Area	Height	Area%
1	7.577	5785	693	0.043
2	8.563	76943	3483	0.577
3	8.748	13124034	1321567	98.355
4	9.046	136807	10196	1.025
Total		13343569	1335939	100.000

邮箱 : marketing@motifbiotech.com

网址 : <https://www.motifbiotech.com>

地址 : 苏州市高新区浒墅关镇浒青路 111 号 B04