


Sample Information

Name : 棕榈酰六肽-12
Sequence : Pal-VGVAPG
Pump A : 0.1% Trifluoroacetic in 100% water
Pump B : 0.1% Trifluoroacetic in 100% acetonitrile
Total Flow : 1ml/min
Wavelength : 214nm
Analytical column type : Phenomenex Luna® 5 µm C18(2) (4.6*250mm*5 µm)
Dissolution method : 15% ACN+85% H₂O
Inj. Volume : 15ul

Time	Pumps	Action	Value
0.01	Module	B.Conc	75
25.00	Pumps	B.Conc	100
33.00	Pumps	B.Conc	100
38.00	Pumps	B.Conc	100
40.00	Pumps	B.Conc	75
50.00	Controller	Stop	

Chromatogram

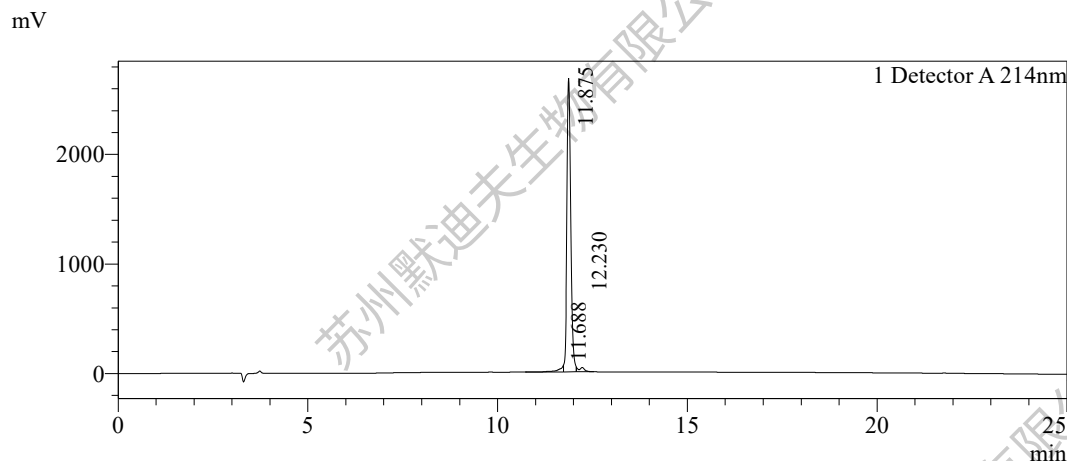

Peak Table

Detector A 214nm

Peak#	Ret. Time	Area	Height	Area%
1	11.688	479905	34180	2.390
2	11.875	19165650	2680725	95.443
3	12.230	435229	40244	2.167
Total		20080785	2755148	100.000

邮箱 : marketing@motifbiotech.com

网址 : <https://www.motifbiotech.com>

地址 : 苏州市高新区浒墅关镇浒青路 111号 B04