


Sample Information

Name : 棕榈酰三肽-1
Sequence : Pal-GHK (HCl Salt)
Pump A : 0.1%Trifluoroacetic in 100% water
Pump B : 0.1%Trifluoroacetic in 100% acetonitrile
Total Flow : 1ml/min
Wavelength : 214nm
Analytical column type : YMC-Triart C18(4.6*250mm*5um)
Dissolution method : 50%ACN+50%H2O
Inj. Volume : 25ul

Time	Module	Action	Value
0.01	Pumps	B.Conc	40
25.00	Pumps	B.Conc	65
33.00	Pumps	B.Conc	100
38.00	Pumps	B.Conc	100
40.00	Pumps	B.Conc	402
50.00	Controller	Stop	

Chromatogram

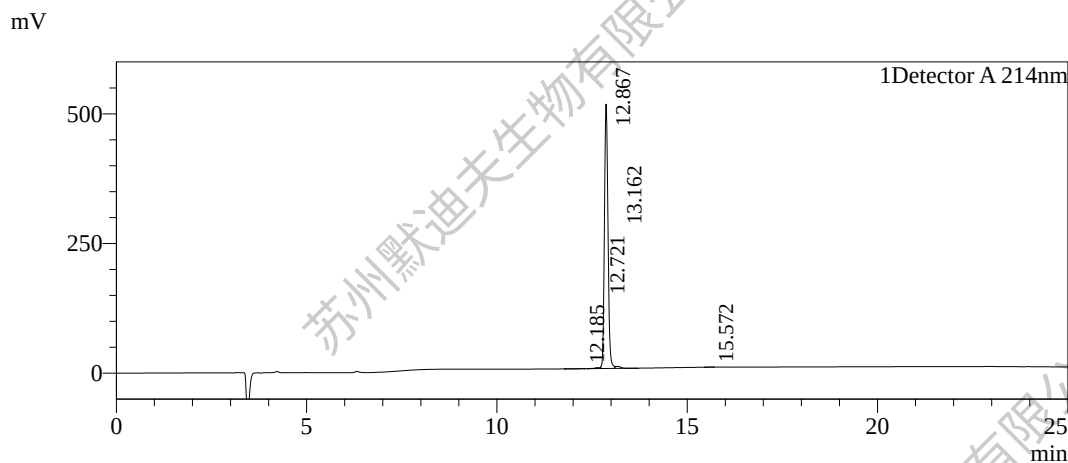

Peak Table

Detector A 214nm

Peak#	Ret. Time	Area	Height	Area%
1	12.185	1711	123	0.051
2	12.721	20479	2225	0.606
3	12.867	3314236	509690	98.141
4	13.162	39002	3813	1.155
5	15.572	1583	233	0.047
Total		3377011	516084	100.000

邮箱 : marketing@motifbiotech.com

网址 : <https://www.motifbiotech.com>

地址 : 苏州市高新区浒墅关镇浒青路 111 号 B04