


Sample Information

Name : 六肽-9
Sequence : GPQGPQ
Pump A : 0.1% Trifluoroacetic in 100% water
Pump B : 0.1% Trifluoroacetic in 100% acetonitrile
Total Flow : 1ml/min
Wavelength : 214nm
Analytical column type : Phenomenex Luna® 5 µm C18(2) (4.6*250mm*5 µm)
Dissolution method : 100% H₂O
Inj. Volume : 25ul

Time	Module	Action	Value
0.01	Pumps	B.Conc	5
25.00	Pumps	B.Conc	30
33.00	Pumps	B.Conc	100
38.00	Pumps	B.Conc	100
40.00	Pumps	B.Conc	5
50.00	Controller	Stop	

Chromatogram

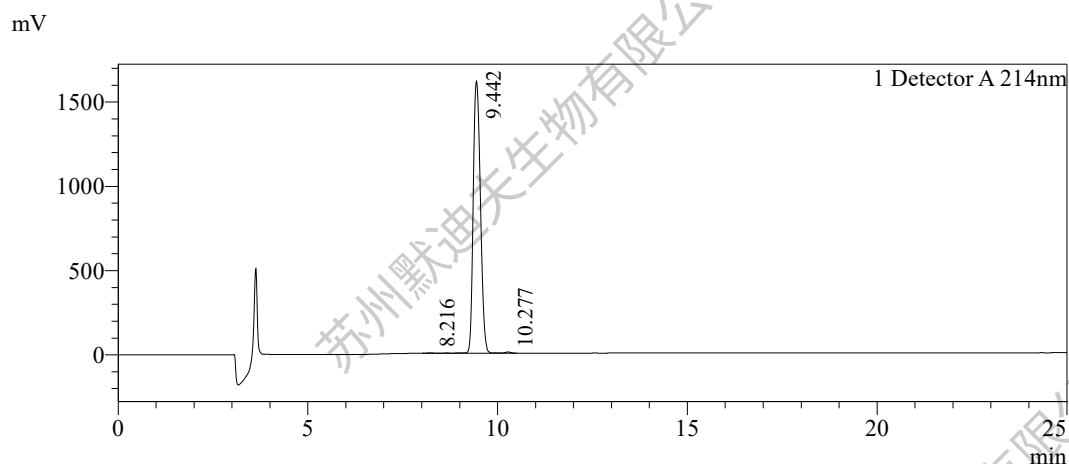

Peak Table

Detector A 214nm

Peak#	Ret. Time	Area	Height	Area%
1	8.216	73059	2584	0.344
2	9.442	21059919	1613735	99.292
3	10.277	77037	6612	0.363
Total		21210015	1622930	100.000

网址 : <https://www.motifbiotech.com>

邮箱 : marketing@motifbiotech.com

地址 : 苏州市高新区浒墅关镇浒青111号 B04